

This is a special file, named RPTHEAD.TXT, in the directory of a method which allows you to customize the report header page.
It can be used to identify the laboratory which uses the method.

This file is printed on the first page with the report styles:

Header+Short, GLP+Short, GLP+Detail, Short+Spec, Detail+Spec, Full

```

      XXXX   XXX
    XX  XX   XX
  XX      XX      XXXXX   XXX XX
  XX      XX XXX   XX    X   XX X XX
  XX    X   XXX XX XXXXXXXX XX X XX
    XX XX   XX XX XX      XX   XX
      XXXX   XXX   XXX   XXXXX   XXX   XXX

```

```

  XXXXXX      X              X      XX
XX   X   XX              XX
XX      XXXXX   XXXXX   XXXXX   XXX   XXXX   XX XXX
  XXXXX   XX      X   XX   XX   XX   XX   XX   XXX XX
    XX   XX   XXXXXX   XX   XX   XX   XX   XX   XX XX
X   XX   XX XX X   XX   XX XX   XX   XX   XX   XX XX
XXXXXX      XXX   XXXXX X   XXX   XXXX   XXXX   XX XX

```

```

                                     X
XX XXX   XXXXX   XX XXX   XXXX   XX XXX   XXXXX
  XXX XX  XX    X   XX XX   XX XX   XXX XX   XX
  XX      XXXXXXXX  XX XX   XX XX   XX      XX
  XX      XX      XXXXX   XX XX   XX      XX XX
XXXXX      XXXXX   XX      XXXX   XXXX      XXX
                XXXX

```

```

  XXX              XXX
  XX              XX
  XX      XXXXX   XXXXX   XX   XXXXX   XX XXX
  XX XXX  XX    X      X   XXXXX   XX    X   XXX XX
  XXX XX  XXXXXXXX  XXXXXXX  XX XX   XXXXXXXX  XX
  XX XX  XX      X   XX  XX XX   XX      XX
  XXX   XXX  XXXXX   XXXXX X   XXXX X   XXXXX   XXXX

```

```

  X              XXX              X
  XX            XX              XX
XXXXX   XXXXX   XXX XX   XX XXX   XX   XXXXX   XXXXX   XXXXX
  XX   XX    X  XX X XX   XX XX   XX      X   XX   XX    X
  XX   XXXXXXXX  XX X XX   XX XX   XX   XXXXXXXX  XX   XXXXXXXX
  XX XX  XX      XX XX   XXXXX   XX   X   XX   XX XX   XX
    XXX   XXXXX   XXX   XXX  XX   XXXX   XXXXX X   XXX   XXXXX
                XXXX

```

```
=====
Acq. Operator   :                               Seq. Line :    5
Acq. Instrument : Kortemmelab HPLC              Location  : Vial 3
Injection Date  : 11/18/2020 12:14:26 AM         Inj       :    1
                                           Inj Volume : 100.0 µl

Sequence File   : C:\Chem32\1\DATA\20201117_GSP1_S75_WITH_STANDAR 2020-11-17 21-00-45\
                  20201117_GSP1_S75_WITH_STANDAR.S
Method          : C:\CHEM32\1\DATA\20201117_GSP1_S75_WITH_STANDAR 2020-11-17 21-00-45\CJM_
                  S75_RUN_ISOCRATIC.M (Sequence Method)
Last changed    : 11/12/2020 5:53:30 PM
Method Info     : S75 analytical
=====
```

```
=====
Module                                     Type   Firmware rev.   Serial number
-----|-----|-----|-----
1200 Autosampler SL                       G1329B A.06.54 [003] DE64155932
1200 Multiple Wavelength Detector          G1365D B.06.72 [0002] DE64256327
1100/1200 Quaternary Pump                 G1311A A.06.32 [011] DE62971812
1200 Sample Thermostat                    G1330B n/a          DEBAK15882
=====
```

Software Revision: Rev. B.04.03 [16] Copyright © Agilent Technologies

=====

Column(s)

=====

```
Column Description : Sephadex 75, 10/300
Serial#            : 10108795
Product#           : 17-5174-01      Batch# :
Diameter           : 4.6 mm          Length : 150.0 mm
Particle size      : 5.0 µm          Void volume : 60.0 %
# Injections       : 75
Maximum Pressure   : 18.0 bar         Maximum pH : 9.0
  Minimum pH      : 2.0
Maximum Temperature: 60.0 °C
Comment           :
=====
```

```
Instrument Conditions :      At Start          At Stop
Pressure             :          0.1              15.6 bar
Flow                 :          0.000            0.800 ml/min
```

Detector Lamp Burn Times: Current On-Time Accumulated On-Time


```
Solvent Description :
PMP1 , Solvent A    : ddH2O
PMP1 , Solvent B    : EtOH
PMP1 , Solvent C    : Tris NaCl
PMP1 , Solvent D    :
```

Run Logbook

18 Nov 20 00:47 AM

Logbook File:C:\Chem32\...ITH_STANDARD 2020-11-17 21-00-45\20201117_WT.D\RUN.LOG

| Module | # Event Message | Time | Date |
|----------|---|----------|----------|
| Method | Method started: line# 5 vial# 3 inj# 1 | 00:14:16 | 11/18/20 |
| Method | Instrument running sample Vial 3 | 00:14:16 | 11/18/20 |
| ALS | Air temperature (tray) = 4.0 °C | 00:17:39 | 11/18/20 |
| PUMP | Pressure = 0.1 bar | 00:17:39 | 11/18/20 |
| PUMP | Flow = 0.000 ml/min | 00:17:39 | 11/18/20 |
| PUMP | Pressure = 15.6 bar | 00:47:40 | 11/18/20 |
| Method | Instrument run completed | 00:47:43 | 11/18/20 |
| Method | Saving Method CJM_S75_RUN_ISOCRATIC.M | 00:47:44 | 11/18/20 |
| Method | Saving Method RUN.M | 00:47:45 | 11/18/20 |
| CP Macro | Analyzing rawdata 20201117_WT.D | 00:47:45 | 11/18/20 |
| CP Macro | Can't load MWD1 D, Sig=230,8 Ref=360,100 | 00:47:45 | 11/18/20 |
| CP Macro | Signal used in Calib. Table (Signal Details)> | 00:47:45 | 11/18/20 |

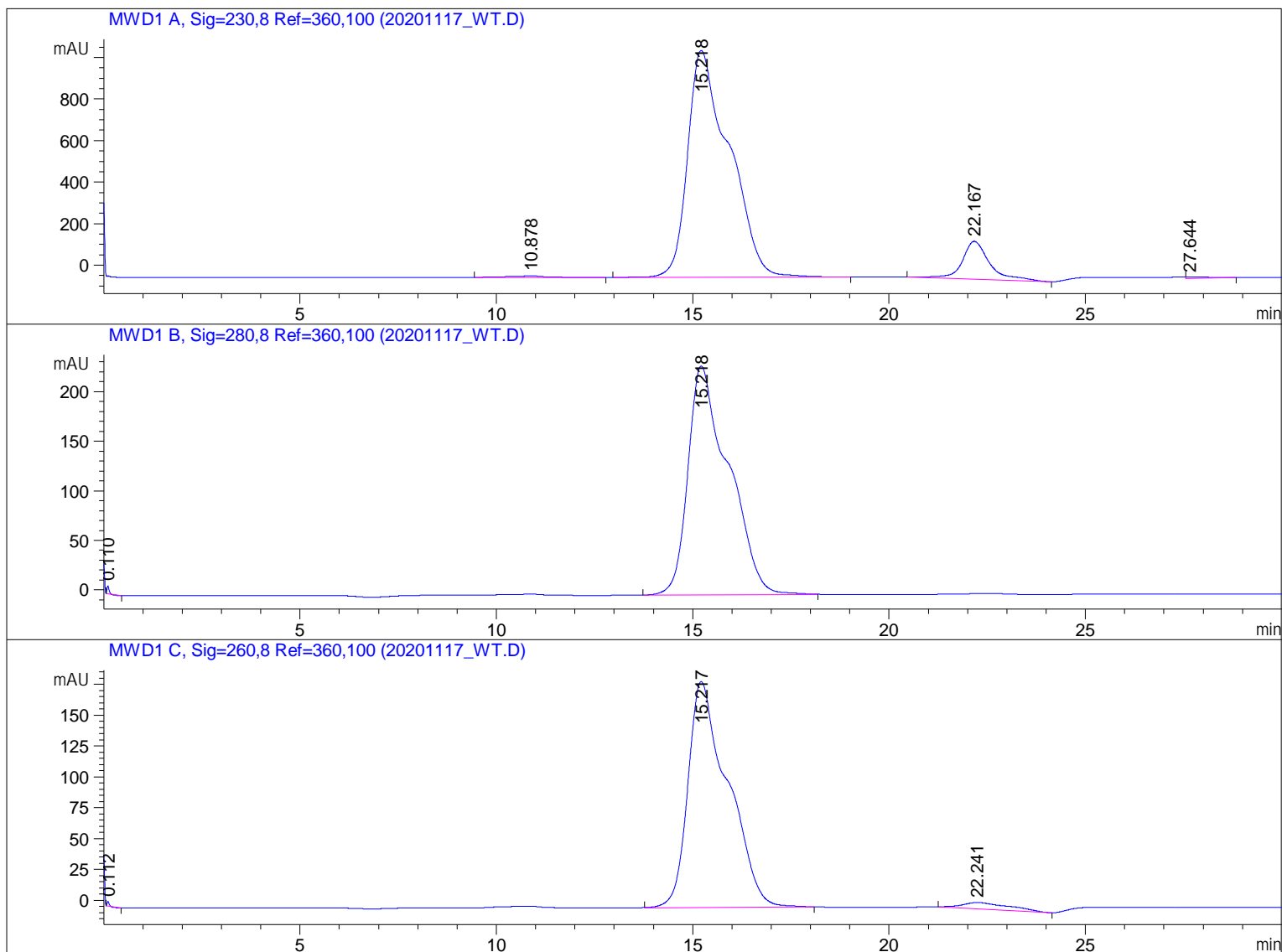

=====
Area Percent Report
=====

Sorted By : Signal
Multiplier: : 1.0000
Dilution: : 1.0000
Use Multiplier & Dilution Factor with ISTDs

Signal 1: MWD1 A, Sig=230,8 Ref=360,100

| Peak # | RetTime [min] | Type | Width [min] | Area [mAU*s] | Height [mAU] | Area % |
|--------|---------------|------|-------------|--------------|--------------|---------|
| 1 | 10.878 | BB | 0.8651 | 506.65463 | 7.72568 | 0.5640 |
| 2 | 15.218 | BB | 1.0241 | 8.00686e4 | 1091.19031 | 89.1262 |
| 3 | 22.167 | BV | 0.7153 | 8943.19727 | 183.26068 | 9.9549 |
| 4 | 27.644 | BB | 0.4908 | 318.90854 | 7.70572 | 0.3550 |

Totals : 8.98374e4 1289.88239

Signal 2: MWD1 B, Sig=280,8 Ref=360,100

| Peak # | RetTime [min] | Type | Width [min] | Area [mAU*s] | Height [mAU] | Area % |
|--------|---------------|------|-------------|--------------|--------------|---------|
| 1 | 0.110 | BB | 0.0420 | 18.90106 | 7.88541 | 0.1143 |
| 2 | 15.218 | BB | 0.9973 | 1.65214e4 | 230.76826 | 99.8857 |

Totals : 1.65403e4 238.65368

Signal 3: MWD1 C, Sig=260,8 Ref=360,100

| Peak # | RetTime [min] | Type | Width [min] | Area [mAU*s] | Height [mAU] | Area % |
|--------|---------------|------|-------------|--------------|--------------|---------|
| 1 | 0.112 | BB | 0.0319 | 6.08425 | 3.56379 | 0.0454 |
| 2 | 15.217 | BB | 0.9853 | 1.29245e4 | 183.14738 | 96.5388 |
| 3 | 22.241 | BV | 1.1290 | 457.29987 | 5.22550 | 3.4158 |

Totals : 1.33879e4 191.93667

=====
*** End of Report ***
=====