

This is a special file, named RPTHEAD.TXT, in the directory of a method which allows you to customize the report header page.
It can be used to identify the laboratory which uses the method.

This file is printed on the first page with the report styles:

Header+Short, GLP+Short, GLP+Detail, Short+Spec, Detail+Spec, Full

```

      XXXX  XXX
    XX  XX  XX
  XX      XX      XXXXX  XXX XX
  XX      XX XXX  XX    X  XX X XX
  XX    X  XXX XX  XXXXXXXX  XX X XX
    XX  XX  XX  XX  XX      XX  XX
      XXXX  XXX  XXX  XXXXX  XXX  XXX

```

```

  XXXXXX      X              X      XX
XX    X  XX              XX
XX      XXXXX  XXXXX  XXXXX  XXX      XXXX  XX XXX
  XXXXX  XX      X  XX      XX      XX  XX  XXX XX
      XX  XX      XXXXXX  XX      XX      XX  XX  XX XX
X    XX  XX XX  X  XX  XX XX      XX      XX  XX  XX XX
XXXXXX      XXX  XXXXX X  XXX  XXXX  XXXX  XX  XX

```

```

                                     X
  XX XXX  XXXXX  XX XXX  XXXX  XX XXX  XXXXX
  XXX XX  XX    X  XX  XX  XX  XX  XXX XX  XX
  XX      XXXXXXXX  XX  XX  XX  XX  XX      XX
  XX      XX      XXXXX  XX  XX  XX      XX XX
XXXXX      XXXXX  XX      XXXX  XXXX      XXX
                XXXX

```

```

  XXX              XXX
  XX              XX
  XX      XXXXX  XXXXX  XX      XXXXX  XX XXX
  XX XXX  XX    X      X  XXXXX  XX    X  XXX XX
  XXX XX  XXXXXXXX  XXXXXXXX  XX  XX  XXXXXXXX  XX
  XX  XX  XX      X  XX  XX  XX  XX      XX
  XXX  XXX  XXXXX  XXXXX X  XXXX X  XXXXX  XXXX

```

```

  X              XXX              X
  XX              XX              XX
XXXXX  XXXXX  XXX XX  XX XXX  XX  XXXXX  XXXXX  XXXXX
XX    XX    X  XX X XX  XX  XX  XX      X  XX    XX    X
XX    XXXXXXXX  XX X XX  XX  XX  XX  XXXXXXXX  XX  XXXXXXXX
XX XX  XX      XX  XX  XXXXX  XX  X  XX  XX XX  XX
  XXX  XXXXX  XXX  XXX  XX      XXXX  XXXXX X  XXX  XXXXX
                XXXX

```

```
=====
Acq. Operator   :                               Seq. Line :    9
Acq. Instrument : Kortemmelab HPLC              Location  : Vial 5
Injection Date  : 11/18/2020 2:29:23 AM          Inj       :    1
                                                Inj Volume : 100.0 µl
Sequence File   : C:\Chem32\1\DATA\20201117_GSP1_S75_WITH_STANDARD 2020-11-17 21-00-45\
                  20201117_GSP1_S75_WITH_STANDARD.S
Method          : C:\CHEM32\1\DATA\20201117_GSP1_S75_WITH_STANDARD 2020-11-17 21-00-45\CJM_
                  S75_RUN_ISOCRATIC.M (Sequence Method)
Last changed    : 11/12/2020 5:53:30 PM
Method Info     : S75 analytical
=====
```

```
=====
Module                                     Type   Firmware rev.   Serial number
-----|-----|-----|-----
1200 Autosampler SL                       G1329B A.06.54 [003] DE64155932
1200 Multiple Wavelength Detector          G1365D B.06.72 [0002] DE64256327
1100/1200 Quaternary Pump                  G1311A A.06.32 [011] DE62971812
1200 Sample Thermostat                     G1330B n/a          DEBAK15882
=====
```

Software Revision: Rev. B.04.03 [16] Copyright © Agilent Technologies

=====

Column(s)

=====

```
Column Description : Sephadex 75, 10/300
Serial#            : 10108795
Product#           : 17-5174-01      Batch# :
Diameter           : 4.6 mm          Length : 150.0 mm
Particle size      : 5.0 µm          Void volume : 60.0 %
# Injections       : 79
Maximum Pressure   : 18.0 bar         Maximum pH : 9.0
Minimum pH         : 2.0
Maximum Temperature: 60.0 °C
Comment           :
=====
```

```
Instrument Conditions :      At Start          At Stop
Pressure              :          0.0            15.5 bar
Flow                  :          0.000          0.800 ml/min
```

Detector Lamp Burn Times: Current On-Time Accumulated On-Time




```
Solvent Description :
PMP1 , Solvent A    : ddH2O
PMP1 , Solvent B    : EtOH
PMP1 , Solvent C    : Tris NaCl
PMP1 , Solvent D    :
```

Run Logbook

18 Nov 20 03:02 AM

Logbook File:C:\Chem32\...H_STANDARD 2020-11-17 21-00-45\20201117_F28Y.D\RUN.LOG

Module	# Event Message	Time	Date
Method	Method started: line# 9 vial# 5 inj# 1	02:29:13	11/18/20
Method	Instrument running sample Vial 5	02:29:13	11/18/20
ALS	Air temperature (tray) = 4.0 °C	02:32:33	11/18/20
PUMP	Pressure = 0.0 bar	02:32:33	11/18/20
PUMP	Flow = 0.000 ml/min	02:32:33	11/18/20
PUMP	Pressure = 15.5 bar	03:02:34	11/18/20
Method	Instrument run completed	03:02:41	11/18/20
Method	Saving Method CJM_S75_RUN_ISOCRATIC.M	03:02:45	11/18/20
Method	Saving Method RUN.M	03:02:47	11/18/20
CP Macro	Analyzing rawdata 20201117_F28Y.D	03:02:47	11/18/20
CP Macro	Can't load MWD1 D, Sig=230,8 Ref=360,100	03:02:47	11/18/20
CP Macro	Signal used in Calib. Table (Signal Details)>	03:02:47	11/18/20

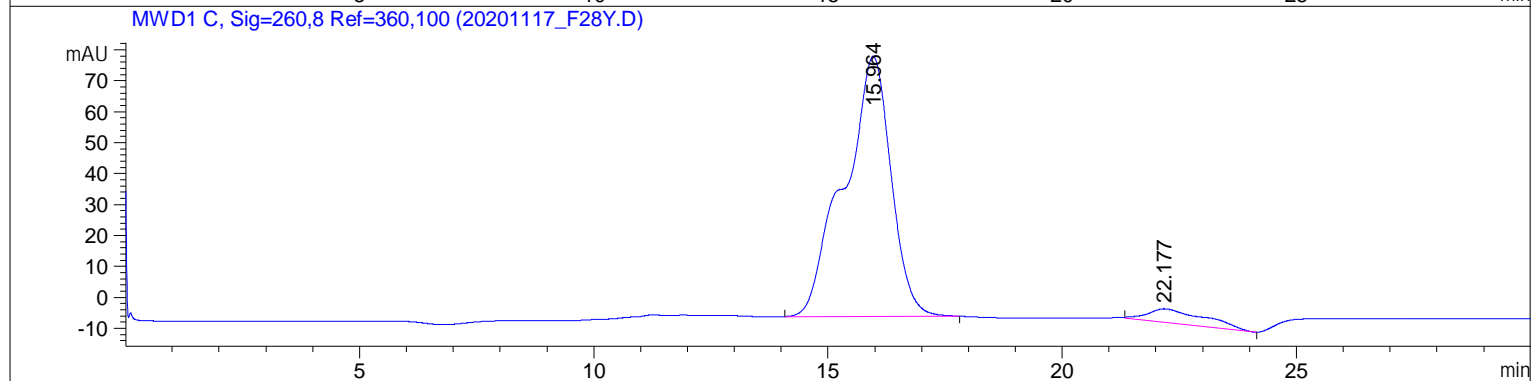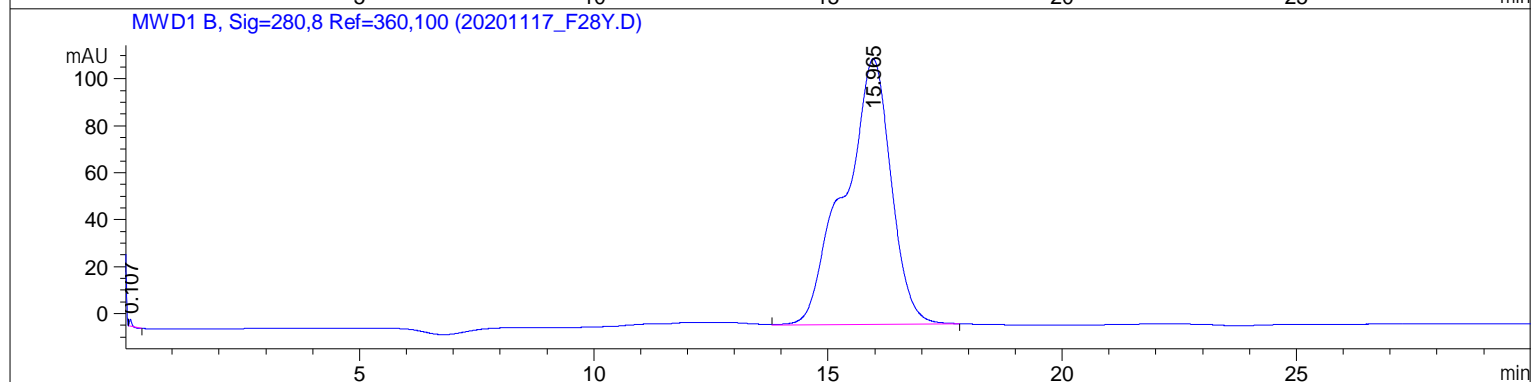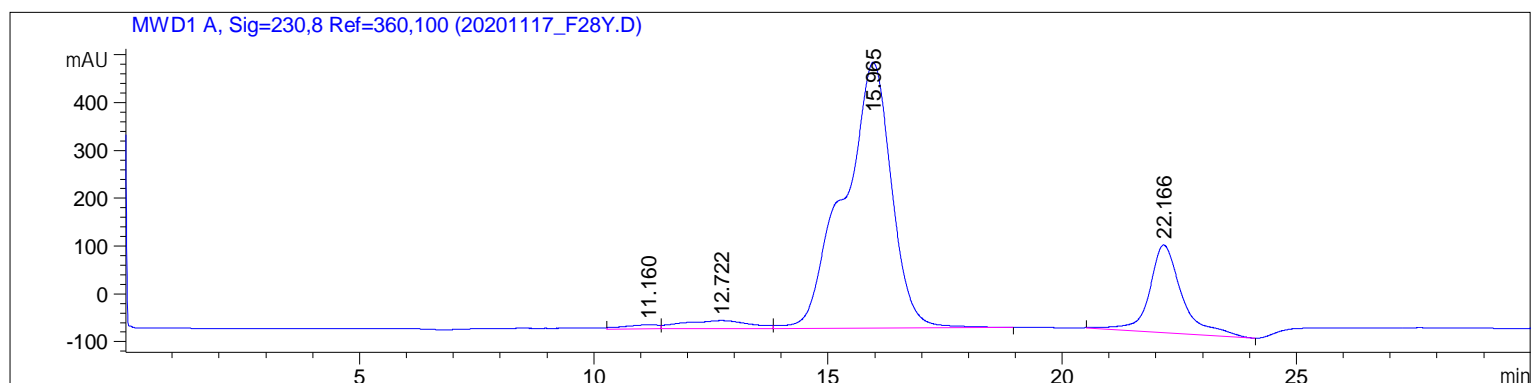

=====
Area Percent Report
=====

Sorted By : Signal
Multiplier: : 1.0000
Dilution: : 1.0000
Use Multiplier & Dilution Factor with ISTDs

Signal 1: MWD1 A, Sig=230,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	11.160	BV	0.6461	434.31744	8.79519	0.8636
2	12.722	VV	1.2542	1679.60168	16.80311	3.3397
3	15.965	VB	1.0016	3.93335e4	555.79749	78.2093
4	22.166	BV	0.7041	8845.22949	183.57545	17.5875

Totals : 5.02927e4 764.97124

Signal 2: MWD1 B, Sig=280,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	0.107	BB	0.0398	6.55403	2.75374	0.0833
2	15.965	VB	0.9816	7857.03271	113.22075	99.9167

Totals : 7863.58675 115.97449

Signal 3: MWD1 C, Sig=260,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	15.964	BB	0.9872	5867.55615	84.17684	93.7981
2	22.177	BV	1.1282	387.95786	4.35117	6.2019

Totals : 6255.51401 88.52802

=====
*** End of Report ***