

This is a special file, named RPTHEAD.TXT, in the directory of a method which allows you to customize the report header page.
It can be used to identify the laboratory which uses the method.

This file is printed on the first page with the report styles:

Header+Short, GLP+Short, GLP+Detail, Short+Spec, Detail+Spec, Full

```

      XXXX  XXX
    XX  XX  XX
  XX      XX      XXXXX  XXX XX
  XX      XX XXX  XX    X  XX X XX
  XX    X  XXX XX  XXXXXXXX  XX X XX
    XX  XX  XX  XX  XX      XX  XX
      XXXX  XXX  XXX  XXXXX  XXX  XXX

```

```

  XXXXXX  X      X      XX
XX  X  XX      XX
XX      XXXXX  XXXXX  XXXXX  XXX  XXXX  XX XXX
  XXXXX  XX      X  XX  XX  XX  XX  XXX XX
    XX  XX  XXXXXX  XX  XX  XX  XX  XX  XX
X  XX  XX XX  X  XX  XX XX  XX  XX  XX  XX
XXXXXX  XXX  XXXXX X  XXX  XXXX  XXXX  XX  XX

```

```

                                     X
  XX XXX  XXXXX  XX XXX  XXXX  XX XXX  XXXXX
  XXX XX  XX  X  XX  XX  XX  XX  XXX XX  XX
  XX      XXXXXXXX  XX  XX  XX  XX  XX
  XX      XX      XXXXX  XX  XX  XX  XX  XX
XXXXX      XXXXX  XX      XXXX  XXXX  XXX
                XXXX

```

```

  XXX      XXX
  XX      XX
  XX      XXXXX  XXXXX  XX  XXXXX  XX XXX
  XX XXX  XX  X  X  XXXXX  XX  X  XXX XX
  XXX XX  XXXXXXXX  XXXXXXXX  XX  XX  XXXXXXXX  XX
  XX  XX  XX      X  XX  XX  XX  XX  XX
  XXX  XXX  XXXXX  XXXXX X  XXXX X  XXXXX  XXXX

```

```

  X      XXX      X
  XX      XX      XX
XXXXXX  XXXXX  XXX XX  XX XXX  XX  XXXXX  XXXXX  XXXXX
  XX  XX  X  XX X XX  XX  XX  XX  X  XX  XX  X
  XX  XXXXXXXX  XX X XX  XX  XX  XX  XXXXXXXX  XX  XXXXXXXX
  XX XX  XX  XX  XX  XXXXX  XX  X  XX  XX XX  XX
    XXX  XXXXX  XXX  XXX  XX  XXXX  XXXXX X  XXX  XXXXX
                XXXX

```

```
=====
Acq. Operator   :                               Seq. Line :    1
Acq. Instrument : Kortemmelab HPLC              Location  : Vial 1
Injection Date  : 11/12/2020 6:26:20 PM          Inj       :    1
                                                Inj Volume : 100.0 µl
Sequence File   : C:\Chem32\1\DATA\20201112_GSP1_F28V_F28Y_S75 2020-11-12 18-26-05\20201112_
                  GSP1_F28V_F28Y_S75.S
Method          : C:\CHEM32\1\DATA\20201112_GSP1_F28V_F28Y_S75 2020-11-12 18-26-05\CJM_S75_
                  RUN_ISOCRATIC.M (Sequence Method)
Last changed    : 11/12/2020 5:53:30 PM
Method Info     : S75 analytical
=====
```

```
=====
Module                                Type   Firmware rev.   Serial number
-----|-----|-----|-----
1200 Autosampler SL                   G1329B A.06.54 [003] DE64155932
1200 Multiple Wavelength Detector      G1365D B.06.72 [0002] DE64256327
1100/1200 Quaternary Pump             G1311A A.06.32 [011] DE62971812
1200 Sample Thermostat                 G1330B n/a          DEBAK15882
=====
```

Software Revision: Rev. B.04.03 [16] Copyright © Agilent Technologies

```
=====
                        Column(s)
=====
```

```
Column Description : Sephadex 75, 10/300
Serial#            : 10108795
Product#           : 17-5174-01      Batch# :
Diameter           : 4.6 mm          Length : 150.0 mm
Particle size      : 5.0 µm          Void volume : 60.0 %
# Injections       : 65
Maximum Pressure   : 18.0 bar         Maximum pH : 9.0
Minimum pH        : 2.0
Maximum Temperature: 60.0 °C
Comment           :
=====
```

```
Instrument Conditions :      At Start          At Stop
Pressure             :          0.0              15.7 bar
Flow                 :          0.000            0.800 ml/min
```

Detector Lamp Burn Times: Current On-Time Accumulated On-Time


```
Solvent Description :
PMP1 , Solvent A    : ddH2O
PMP1 , Solvent B    : EtOH
PMP1 , Solvent C    : Tris NaCl
PMP1 , Solvent D    :
```

Run Logbook

12 Nov 20 06:58 PM

Logbook File:C:\Chem32\...28Y_S75 2020-11-12 18-26-05\20201112_BLANK_1.D\RUN.LOG

Module	# Event Message	Time	Date
Method	Method started: line# 1 vial# 1 inj# 1	18:26:10	11/12/20
Method	Instrument running sample Vial 1	18:26:10	11/12/20
ALS	Air temperature (tray) = 4.0 °C	18:28:37	11/12/20
PUMP	Pressure = 0.0 bar	18:28:37	11/12/20
PUMP	Flow = 0.000 ml/min	18:28:37	11/12/20
PUMP	Pressure = 15.7 bar	18:58:38	11/12/20
Method	Instrument run completed	18:58:41	11/12/20
Method	Saving Method CJM_S75_RUN_ISOCRATIC.M	18:58:42	11/12/20
Method	Saving Method RUN.M	18:58:44	11/12/20
CP Macro	Analyzing rawdata 20201112_BLANK_1.D	18:58:44	11/12/20
CP Macro	Can't load MWD1 D, Sig=230,8 Ref=360,100	18:58:44	11/12/20
CP Macro	Signal used in Calib. Table (Signal Details)>	18:58:44	11/12/20

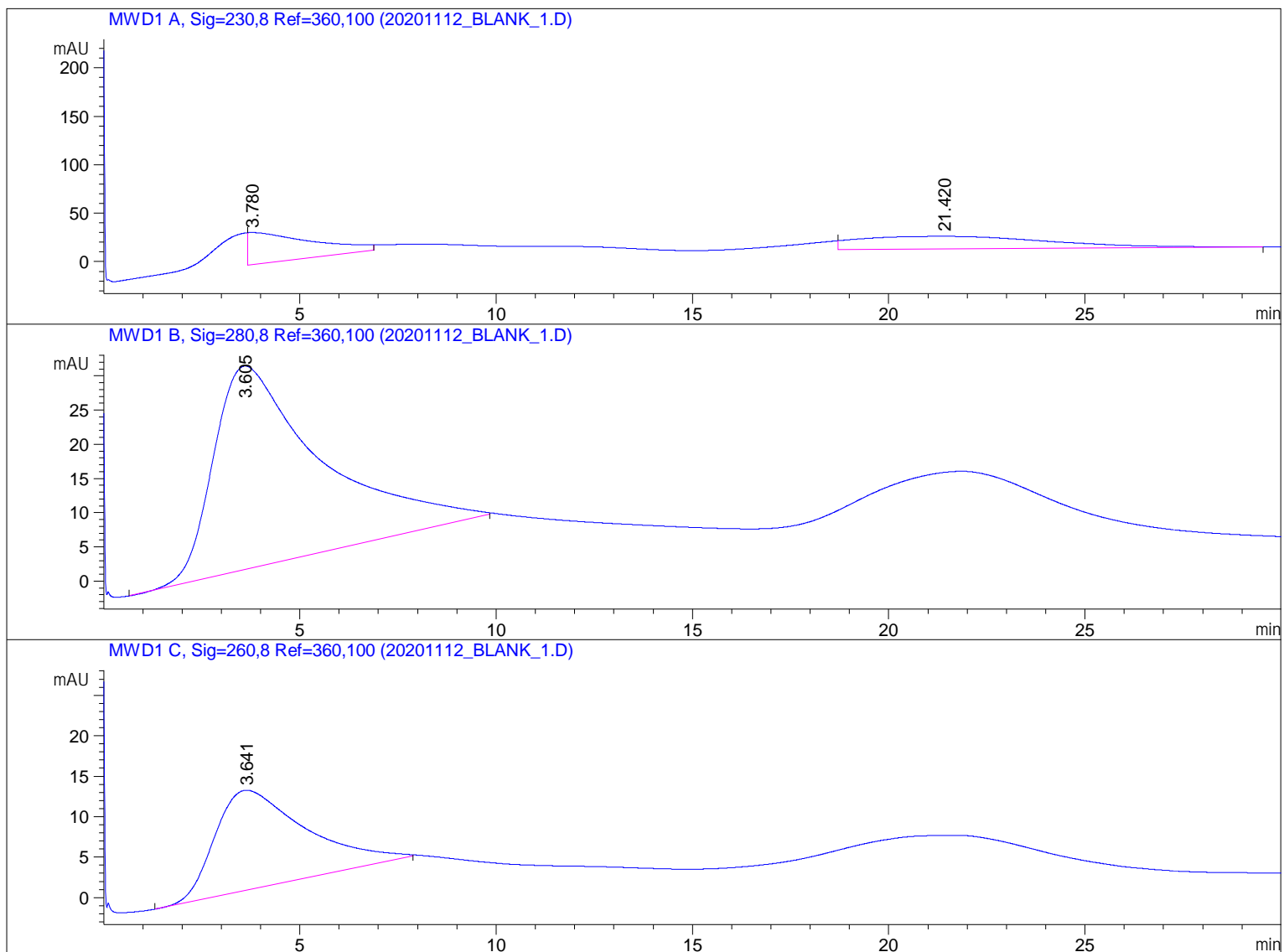

=====
Area Percent Report
=====

Sorted By : Signal
Multiplier: : 1.0000
Dilution: : 1.0000
Use Multiplier & Dilution Factor with ISTDs

Signal 1: MWD1 A, Sig=230,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	3.780	BV	1.2361	3497.74512	33.28162	43.9840
2	21.420	BB	3.9823	4454.56494	13.06996	56.0160

Totals : 7952.31006 46.35159

Signal 2: MWD1 B, Sig=280,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	3.605	BB	2.4646	5623.23438	29.75784	100.0000

Totals : 5623.23438 29.75784

Signal 3: MWD1 C, Sig=260,8 Ref=360,100

Peak #	RetTime [min]	Type	Width [min]	Area [mAU*s]	Height [mAU]	Area %
1	3.641	BB	1.8676	1964.13025	12.36241	100.0000

Totals : 1964.13025 12.36241

=====
*** End of Report ***